


Experiment – 1

Aim: To determine the strength of ammonia in tap water sample.

Apparatus required: Pipette, burette, funnel, conical flask, volumetric flask, and dropper

Chemical required: Standard N/20 - Na_2CO_3 , HCl solution, water sample and methyl orange indicator.

Theory: NH_3 is highly soluble in H_2O . It is nitrogenous compound and present in water in form of ammonia hydroxide in solution.

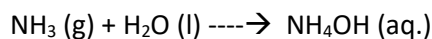

Where, NH_3 is coming from the pesticides, insecticides fertilizer and other industrial water that are disposed off into the water.

The titration of NH_3 and HCl is called Acid – Base titration. Here, HCl is not a primary standard solution, So it is firstly standardized by the help of Na_2CO_3 solution.

Procedure:

I) Standardisation of HCl Solution -

- i- Rinse and fill the burette with HCl.
- ii- Pipette out 10ml of Na_2CO_3 solution and added 2 drops of methyl orange, titrated with HCl solution till colour changed from yellow to light pink.
- iii- The experiment was repeated three times to get concordant reading.

II) Titration of water sample with HCl Solution –

- i- Rinse and fill the burette with HCl.
- ii- Pipette out 10ml water sample in conical flask.
- iii- Add 2-3 drops of methyl orange and titrate with HCl solution till the colour change from yellow to light pink.
- iv- The experiment was repeated three times to get concordant reading.

Observation Table:

I) Standardisation of HCl-

S. No.	Burette Reading		Volume of HCl used, (ml)
	Initial	Final	
1	0	13	13
2			
3			

Concordant reading = ml

$$\begin{aligned} N_1 V_1 &= N_2 V_2 \\ (\text{Na}_2\text{CO}_3) & \quad (\text{HCl}) \\ (1/20) \times 10 &= N_2 \times \text{.....} \\ N_2 &= \text{..... N of HCl} \end{aligned}$$

(II) Titration of water sample with HCl

S. No.	Burette Reading		Volume of HCl used(ml)
	Initial Volume (ml)	Final Volume (ml)	
1	0	2.2	2.2
2			
3			

Concordant reading 2.2 ml

$$\begin{aligned} N_2 V_2 &= N_3 V_3 \\ (\text{HCl}) & \quad (\text{Water}) \\ \dots \times \dots &= N_3 \times 10 \\ N_3 &= \dots \text{ N of water} \end{aligned}$$


Strength = $N_3 \times 17$

$$= \dots \text{ g/L of NH}_3$$

Structure:

Indicator -Methyl Orange

(Acidic Medium) - Pink

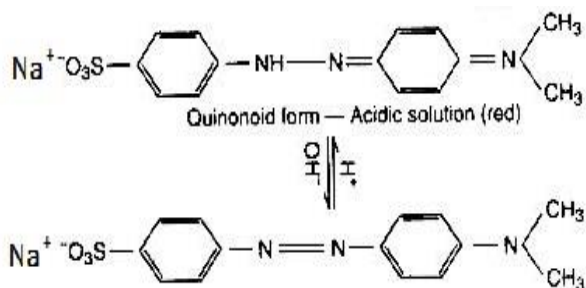

(Basic Medium) - Yellow

Result:

Normality of water =N

Normality of HCl =N

Strength of NH_3 = g/L

Precautions:

1. Handle the apparatus with carefully.
2. Do not pipette out hot solutions.
3. There should be no air bubbles in burette.

4. Lower meniscus of HCl should be observed.